Bilancia pesapersona  
Bathroom Scale  
Pèse-personne  
Personenwaage  
Báscula

**Slimmy 2500 Medical**

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**ENGLISH****CORRECT PRODUCT USE**

This appliance must be used only for the use for which it is intended: as a bathroom scale.

Any other use is incorrect and therefore hazardous. Polti S.p.A. may not be held liable for any damage or injury resulting from incorrect or mistaken use.

**WARNINGS**

Read the instructions carefully before using the appliance and keep them for future reference.

This appliance is intended exclusively for domestic use, and must not be used in hospitals or healthcare facilities.

If the device is faulty or malfunctioning, do not attempt to dismantle it yourself, but contact an authorised Polti service centre.

If the appliance is dropped, have it inspected at an authorised service centre.

Wipe with a damp cloth after use. Do not use chemical products or submerge it in water.

Do not subject the appliance to excessive vibration or blows. Do not drop it onto the floor.

Do not leave the appliance where it will be exposed to the weather.

Do not leave the appliance near a heat source.

**Features**

- Capacity 160 Kg.

**Operation**

Place the scale on a flat, solid surface (avoid carpets).

Adjust the scale before use with the wheel on the back of the scale.

Stand still on the scale while you are being weighed, staying in the centre of the scale and waiting for the scale to calculate your weight.


**MAINTENANCE**

Clean the surface of the scale with a damp cloth, but make absolutely sure no water gets into the scale. Do not use chemical products or other cleaning products.

Treat the scale carefully as it is a high precision instrument. Do not drop it or jump on it.

**IMPORTANT**

- This appliance is manufactured in compliance with current recycling regulations.
- Do not dispose of protective packaging material with household waste, but take it to a collection facility.
- Use only original Polti spare parts.

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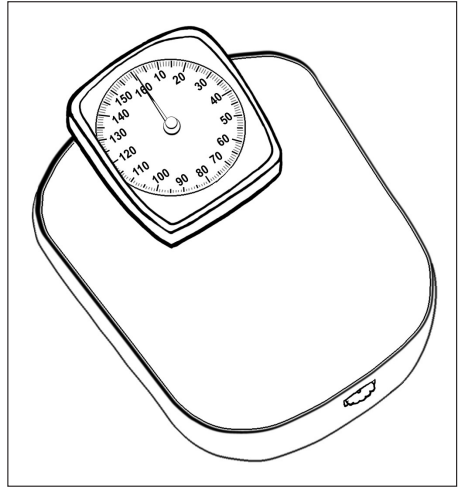**ENGLISH**

Polti products are guaranteed to the original purchaser for two years from the date of purchase as shown on a valid receipt issued by the seller.

This guarantee does not affect the consumer's rights under European Directive 99/44/CE regarding certain aspects of the sale and guarantee of consumer goods; consumers should address claims under these provisions to the retailer where they purchased the goods.

This guarantee is effective in countries where European Directive 99/44/CE has been implemented. Local guarantee legislation applies in all other countries.

During the guarantee period, Polti will repair products with defects in workmanship or materials at no charge to the customer for either labour or materials.

To have work performed under this guarantee, consumers must contact an Authorised Polti Technical Assistance Centre. A list of these

centres and their addresses is attached. The consumer must present a valid receipt issued by the seller at the time of purchase as evidence of the date of purchase.

The guarantee does not cover defects due to incorrect use or use other than that identified in the instruction booklet, which is an integral part of the contract of sale for the product. The guarantee does not cover defects resulting from chance (fire, short circuit) or attributable to a third party (tampering with the product).

The guarantee does not cover parts that are damaged by wear or subject to normal wear (consumables such as filters, brushes, etc.).

Polti will accept no responsibility for any injury to people or animals or damage to objects directly or indirectly resulting from failure to follow the instructions provided in the instruction booklet, in particular the warnings and product maintenance instructions.